



BOWIE&DICK TEST PACK

DAILY AIR REMOVAL TEST [Code: BD3090]




Before



After



Front

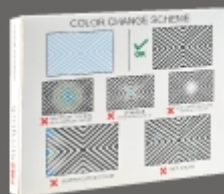

Back

Air leaks in the sterilization chamber threatens the sterilization safety. Daily use of Bowie&Dick test packs helps to determine sterilizer malfunctions at an early stage.

The reactive ink gives a warning before the load leaves CSSD. PMSSteripack Bowie&Dick Test packs comes preassembled, ready to use. No preparation is needed.

The process indicator on the packs' label shows the pack has been used. The test cards can be stored for record purpose.



STEAM INDICATOR STRIP

[Code: IS2001]



- Dual strip can be divided into two for economy of use
- Conforms with ISO 11140-1
- Changes color to BLACK after successful Steam sterilization cycle
- Supplied in boxes of 250 dual



STEAM INDICATOR STRIP

[Code: IS2002]



- Conforms with ISO 11140-1
- Changes color to GREEN after successful Steam sterilization cycle
- Supplied in boxes of 250 pcs.



EO INDICATOR STRIP

[Code: IE3001]



- Dual strip can be divided into two for economy of use
- Conforms with ISO 11140-1
- Changes color to YELLOW after successful EO sterilization cycle
- Supplied in boxes of 125 dual



VH2O2 INDICATOR STRIP

[Code: IP4001]



- Conforms with ISO 11140-1
- Changes color to YELLOW/ORANGE after successful sterilization cycle
- Supplied in boxes of 50 pcs.